


MICROFIL[®] Injection Compounds

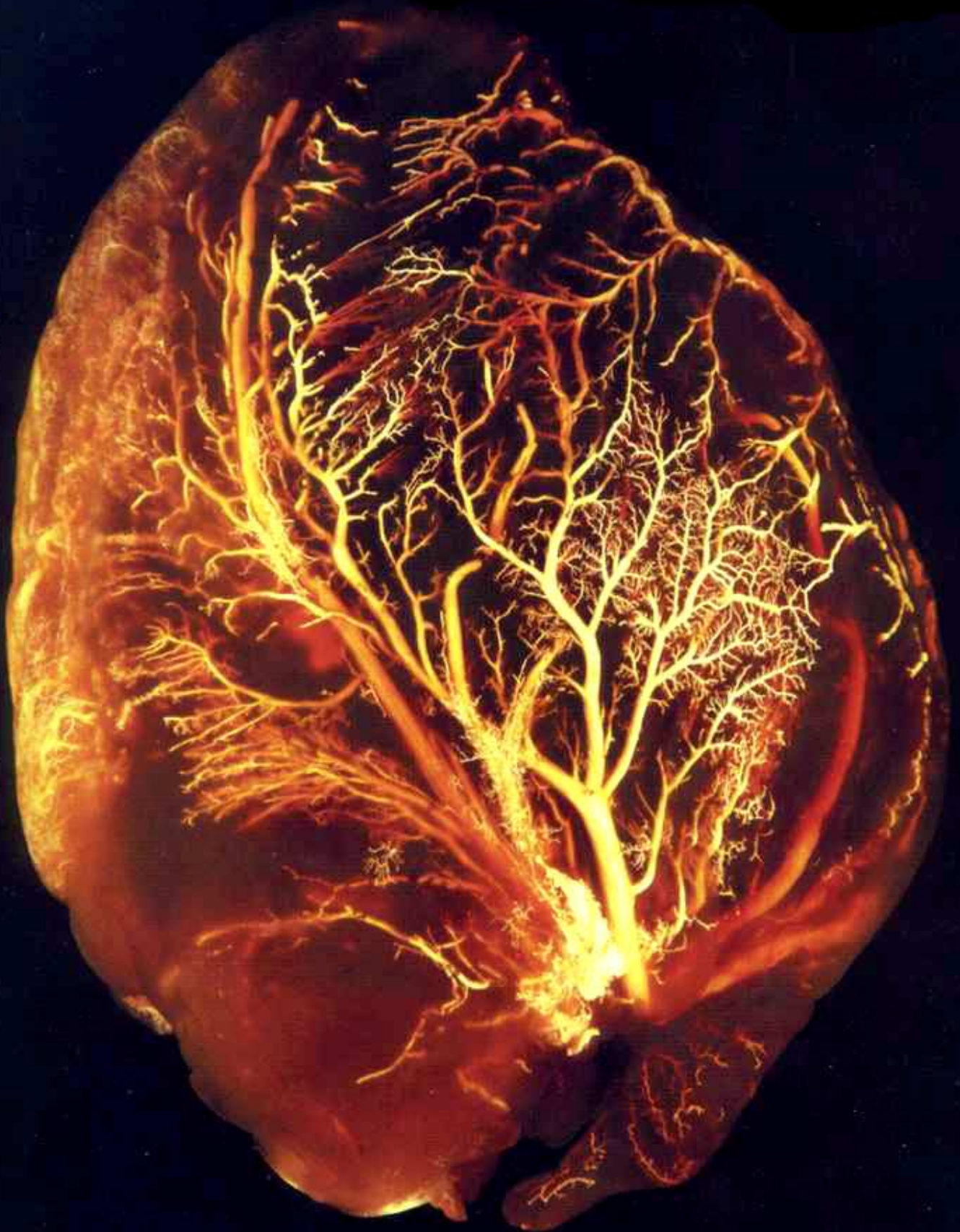

MICROFIL® compounds will fill and opacify microvascular and other spaces of non-surviving animals and postmortem tissue under physiological injection pressure. The continuous, closed vascular system tends itself to flow through injection or perfusion techniques. Following injection, MICROFIL compounds cure to form a three-dimensional cast of the vasculature.

MICROFIL MV-series compounds are available in five radiopaque colors, as well as clear. MV-series compounds require either an alcohol-methyl salicylate or glycerin clearing sequence, whereby the refractive index of the clearing solution is the same as the refractive index of the tissue. This allows for microscopic examination of a selected vascular bed.

MICROFIL CP-101 compound is intended for use in cast-corrosion techniques, and is designed for filling large blood vessels (greater than 100 microns). Although the CP-101 compound will fill capillaries, these vessels fragment when the supporting tissue is removed through exposure to a potassium hydroxide solution. When cured, MICROFIL CP-101 is milky white in color. Casts made using CP-101 will maintain their dimensional accuracy indefinitely.

Advantages offered with MICROFIL compounds over previously available rubber injection materials include:

- Complete filling with minimal shrinkage, to enhance vessel continuity and to produce in cleared preparations a vivid, optically cleared specimen that allows a precise study of the microcirculation.

- Color diversity, to provide delineation within the circulatory tree for microscopic examination and photographic illustration.

Areas of investigation

In physiology, MICROFIL visualization provides a means for establishing the precise vascular architecture of specific organs, allowing comparison between normal and abnormal structure.

In surgery, visualization of the microcirculation and microanatomy is leading to improved surgical techniques in the repair of nerves, tendons, and blood vessels.







In gastrointestinal research, MICROFIL compounds characterize and describe changes in vascular patterns associated with several pathological conditions.

Injected specimens, when preserved in methyl salicylate or glycerin, also serve as a definitive teaching adjunct.

MV-series mixing procedure

To achieve a viscosity level suitable for injection of the microcirculation, it is necessary to blend the MV compound with an equal quantity (by weight) of MV-Diluent. Volume mixing requires 5ml of diluent for every 4ml of compound. The mixture of compound and diluent is catalyzed with 5% (by weight or volume) of MV Curing Agent; Viscosity ranges from 20 to 30 centipoise. Working time is 20 minutes and begins with the addition of curing agent.

Table 1.
Physical properties
of MICROFIL
compounds

	MV-112	MV-117	MV-120	MV-122	NW-130	MV-132	MV-Diluent	CP-101
Color	 White	 Orange	 Blue	 Yellow	 Red	Clear	Clear	 Milky White
Specific gravity	1.04	1.02	1.00	1.04	1.02	1.00	0.92	0.98
Viscosity, ¹ centipoise	35	25	25	25	30	20	5	25
Gel time, ² minutes	90	90	90	90	90	90	—	45
Useful Shelf life, ³ months	4	4	4	4	4	4	Indefinite	6

Notes

1. Viscosity measured with a Brookfield Model LVF Viscometer, with a No. 2 spindle at 30 RPM.
2. Gel time measured on a blend (by weight) of one part MV compound and one part MV-Diluent followed by addition of 5% MV Curing Agent. Gel time is time required for mixture to cease flowing.
3. Shelf life is in excess of four months. For materials held beyond four months, it is possible to run the following static test to determine acceptability:
 - a. Mix 5 grams of MV compound with 5 grams of MV-Diluent in a vial.
 - b. Add 10 drops (medicine dropper) of MV Curing Agent.
 - c. Cap, shake, and refrigerate overnight.
 - d. The mixture should gel in the vial after this procedure to assure a cure in subsequent animal injections.

Catalyzed mixtures will form an elastomeric gel after 90 minutes at room temperature. Curing takes place with non-exothermic cross-linking and minimal volume change.

It has been possible to refrigerate specimens immediately after injection and still obtain complete cure after overnight aging. This procedure decreases odor level for subsequent sectioning.

Perfusion techniques

Two techniques of tissue clearing are described below. Alcohol-methyl salicylate clearing produces a stiffer tissue which, from an aesthetic point, provides a pleasing view for gross observation. Glycerin clearing produces a more flexible tissue, allowing easier manipulation for a given vessel.

Alcohol-methyl salicylate clearing

Non-wetting features of MICROFIL compounds prevent any interaction with blood. Therefore, in the non-surviving animal, a selected vascular bed can be readily perfused without prior washout of blood. Heparinization to maintain blood fluidity, however, has been used to realize improved injection preparations.

For the injection of blood vessels from vascular beds removed postmortem, washout of clotted blood with saline is advisable.

Selected vascular beds are perfused through their accessible artery and drained through a similar vein. Infusion pressures will vary with the animal's mean systemic pressure. For organs from the dog, cat, rat,



Ⓢ MICROFIL MV-130
(Red) injection of rat
trachea.
*Courtesy Robert A. Acland,
University of Louisville*



Ⓢ MICROFIL MV-130
(Red) injection of rat
kidney.
*Courtesy Robert A. Acland,
University of Louisville*

Table 2
Alcohol-methyl salicylate clearing sequence

First Day	—	Immerse in a 25% solution of ethyl alcohol.
Second Day	—	Immerse in a fresh solution of 50% ethyl alcohol.
Third Day	—	Immerse in a fresh solution of 75% ethyl alcohol.
Fourth Day	—	Immerse in a fresh solution of 95% ethyl alcohol.
Fifth Day	—	Immerse in absolute ethyl alcohol.
Sixth Day	—	Immerse for 12 to 24 hours in methyl salicylate.

If tissue has not cleared, return to 95% ethyl alcohol stage and repeat final steps of clearing procedure.

Notes

1. At the 50% ethyl alcohol concentration, tissue specimens may be bleached with 6% hydrogen peroxide for one day. After bleaching, continue the normal clearing procedure. Peroxide bleaching permits greater depth perception; however, this procedure must be considered against some loss in color contrast.
2. In the final stages of ethyl alcohol-methyl salicylate clearing, the clearing liquid may have a cloudy appearance. Addition of a small amount of 10% ethyl alcohol will alleviate this condition.

and man, a pressure of 150mm. Hg for arterial fillings has been used, and 25 to 50mm.Hg for venous fillings.

After the vascular bed is perfused and the MICROFIL injection mass allowed to cure overnight at room temperature, the tissue is subjected to the clearing sequence described in Table 2. Thin tissues may be cleared without sectioning, but thicker organs, such as kidney or brain, should be cut into 1-centimeter slices before immersion. The alcohol-methyl salicylate clearing technique is applicable to all types of tissue with the exception of brain tissue. In the case of brain tissues, it is necessary to allow two days for each step, with an alcohol solution change every day.

Glycerine clearing

The animal is anesthetized with Nembutal, 25 mg/kg i.v., at the same time it is heparinized, to ensure effective removal of its blood volume during perfusion.

A midline incision exposes the abdominal viscera from the sternal notch to symphysis. Following the placement of an abdominal retractor, the thoracic cage is opened rapidly, the thoracic aorta isolated, and a polyethylene cannula inserted distally. The arterial cannula is connected to a sine wave perfusion pump and, prior to instigating perfusion, the right atrium is opened to serve as a drain vent.


The animal is perfused with saline until all of the visceral blood volume is flushed out and the perfusate drained through the arterial vent is essentially free of blood. Adequate perfusion is characterized by severe blanching of all visceral organs.

During perfusion, the curing agent is added to the MICROFIL injection mass. When perfusion is complete, the silicone rubber (e.g., MV-130 Red) is infused through the aortic cannula by syringe. When filling is complete, all organs have a rich, red coloration. MICROFIL compound infusion is continued until the injection mass flows freely from the atrial vent. The atrium and arterial cannula are then clamped and the animal is placed under refrigeration at 4°C overnight, to allow polymerization.

On the following day, specimens are taken by careful dissection, and placed in a 50% mixture of water and glycerin. At successive 24-hour intervals, the glycerin concentration is raised to 75%, then 85%, and finally pure glycerin. This procedure clears the tissue so that microscopic examination readily allows three-dimensional visualization of the vascular bed.

Cross-sectional view of monkey jejunum perfused with MICROFIL MV-118 (Maroon)* followed by MV-122 (Yellow).

*Courtesy D. G. Reynolds,
Walter Reed Army Institute
of Research*

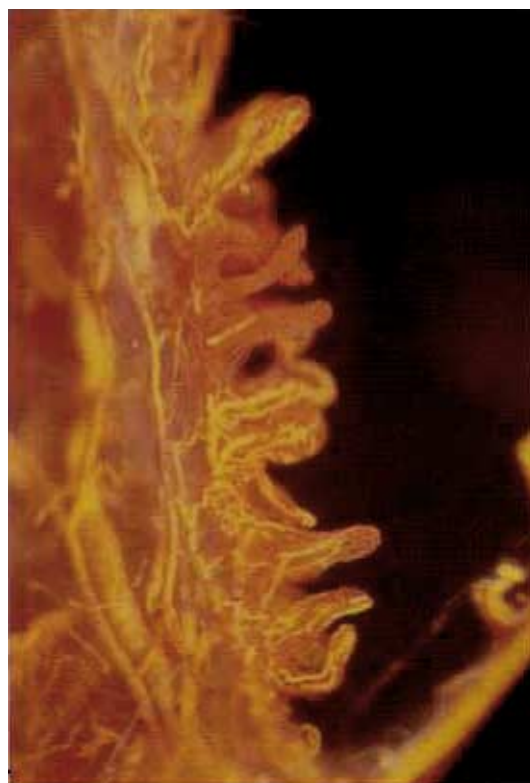

* MV-118 (Maroon) is no longer available and has been replaced with MV-130 (Red).

Additional Notes

1. All MV-series compounds are compatible with one another. Therefore, it is possible to mix colors together to suit your needs (e.g., mix MV-120 Blue with MV-122 Yellow to produce a green compound).
2. Faster cure rates are possible by replacing the conventional MV Curing Agent with 2% ethyl silicate and 1% stannous octoate. Using this cross-link and curing agent combination decreases working time to 5 minutes with complete cure in 20 minutes. If you desire this type of cure system, please specify when placing your order. There is no extra charge for substituting this cure system in place of MV Curing Agent.
3. Occasionally it may be necessary to decrease viscosity by changing the mix ratio to either 2 or 3 parts MV-Diluent for each part MV compound. If your study requires such action, the correct level of MV Curing Agent is 10% (by weight) of the amount of MV compound used.
4. One procedure for vessel differentiation is to completely fill a given circulation with conventional MV compound (e.g., MV-130 Red) then, once the circulation is filled, immediately re-inject the artery with a high-viscosity version of MV compound in a different color (e.g., MV-120 Blue). This is accomplished by substituting MV-Diluent, which has a viscosity of 5 cps, with HV-Diluent, which has a viscosity of 1000 cps. With HV-Diluent, the viscosity of MV-130 Red increased to 350-450 cps. Although this procedure stopped penetration at the capillary level, complete success was obscured by shunting

activity. If you desire to replace MV-Diluent with HV-Diluent, please specify when placing your order. There is no extra charge for this substitution.


5. Store MICROFIL kits at room temperature for maximum retention of pigment dispersion. Keep all containers tightly capped. If pigment settling occurs, it is better to lightly shake the MICROFIL compound container and decant this portion. Stirring the container may put an agglomerated pigment particle into the system, which can be detrimental to perfusion.

6. If catalyzed material should spill on clothing, the best available solvent is MV-Diluent. To facilitate removal of cured compound, swelling and softening will occur on contact with an aromatic solvent such as toluene or xylene, and chlorinated solvents such as trichloroethylene.

Arterial injection of a dog kidney using MICROFIL MV-112 (White). Demonstrates filling of glomeruli and peritubular capillaries at 50x magnification.

Courtesy A. C. Berger, Harvard Medical School

≈

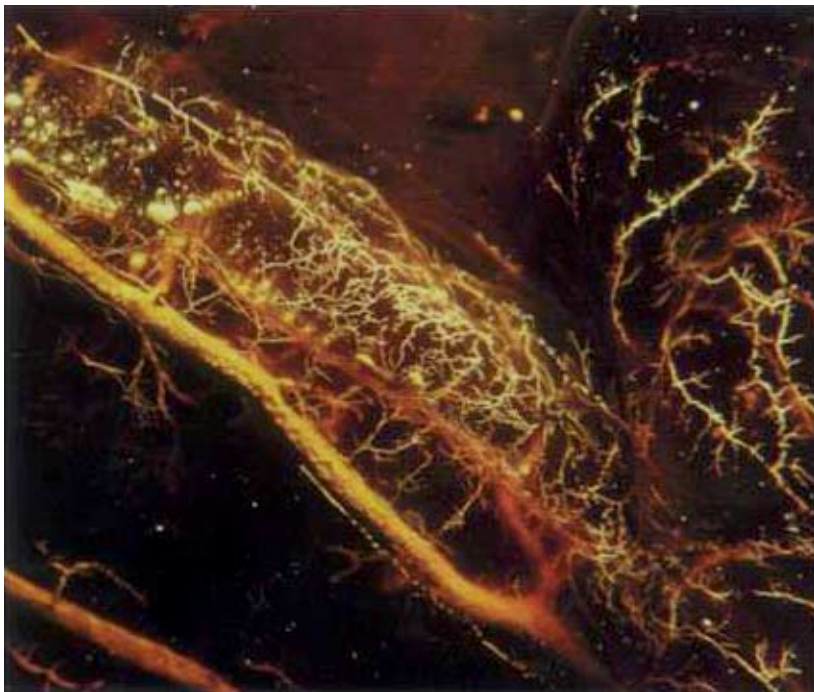

5 Vasa vasorum of the human coronary artery filled with MICROFIL compound.

Courtesy of A. C. Barger, R. Beeuwkes 111, L. L. Lainey and K. J. Silverman, Harvard Medical School

Cover Photo:
MICROFIL injection of the left coronary artery of
the human heart.
*Courtesy of A.C. Barger, R. Beeuwkes III, L.L. Lainey and
K.J. Silverman, Harvard Medical School*

Ordering Information

When ordering material, please remember kit weight is based upon the combined weight of MV compound and MV-Diluent (i.e., a 1-pound kit contains 8 ounces of MV compound and 8 ounces of MV-Diluent). In addition, each kit contains enough MV Curing Agent to cure the contents.

Kit Size Specifications	1 lb.	Any one color
	2 lb.	Any one or two colors
	8 lb.	(1 gallon) Any combination of colors

MICROFIL injection kits are available from:

Flow Tech, Inc.
P.O. Box 834
Carver, Massachusetts 02330
Tel: (508) 866-0007
Fax: (508) 866-0090

Terms: Net 30
FOB Carver

A bibliography is available on request.

While the information herein is believed to be reliable, Flow Tech does not guarantee its accuracy. Users are urged to perform their own tests. MICROFIL products are sold without warranty, and various patents may be pertinent to their use and to the use of compositions containing them. The information contained herein is not intended as a recommendation to use MICROFIL products so as to infringe on any patent. Flow Tech assumes no liability for the user's violation of patent or other rights.

©1996, Flow Tech, Inc. All rights reserved.

MICROFIL is a registered trademark of Flow Tech, Inc., Carver, Massachusetts.